



Surgical perfection. For life.

Setting new low-dose standards in the OR

The STILLE imagiQ2™ surgical imaging table represents an important advance in developing the low-dose technology field, making the OR an even safer place, for patients and for the entire surgical team.


More efficient and shorter procedure times

The True Free Float® technology* of the STILLE imagiQ2 gives unprecedented control and freedom of movement in any desired direction while providing state-of-the-art smoothness like an extension of the surgeon's hand, allowing for shorter procedure times and higher efficiency in the OR.

Choose a low-dose enabler

The unique high-performance carbon fiber table top*, provides unbeatable transparency of 0.4 mm Al equivalency, enabling the C-arm to reduce the radiation dose without compromising image quality. That is 60% more translucency than average tables on the market.

*Patent pending

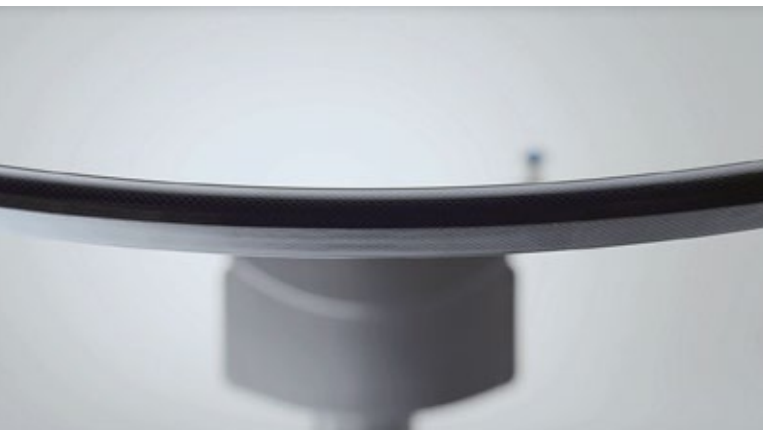

0.4 mm Aluminum Equivalency table top enhancing image quality while reducing radiation exposure for the patient and the entire surgical team.

STILLE imagiQ2 WINNER of PEOPLE'S CHOICE AWARDS 2014



Optimizing C-arm capabilities

With its ultra-slim design and super-low base, the imagiQ2 enhances C-arm access, making it the ideal C-arm complement for the mobile Hybrid OR or office suite. Furthermore, the low-dose table top of the STILLE imagiQ2 contributes to avoiding C-arm overheating.


Continuous image focus

The patented technology of the iso-roll cradling motion ($\pm 15^\circ$) helps maintain image focus even while in full lateral roll, minimizing patient and C-arm repositioning, as well as reducing scattered radiation.

STILLE imagiQ2

- Translucent length: patients as tall as 6'5" / 196 cm can be fully imaged with headrest (5'6" / 169.2 cm without headrest). 360° unobstructed translucent length: 3'6" / 107.2 cm with headrest (2'7" / 80.4 cm without headrest).
- Patient weight capacity: 225 Kg (496 lbs). Can be upgraded to handle patients as heavy as 300 Kg (661lbs) with the *Load Upgrade Pack*.
- True Free Float® allows smooth longitudinal (35.4" / 90 cm), lateral (9.9" / 25 cm) and diagonal (36.8" / 93.5 cm) movements.
- Upgrades to fit your needs: Load upgrade pack, performance pack, mobility enhancement pack.

Learn more at: www.stille.se

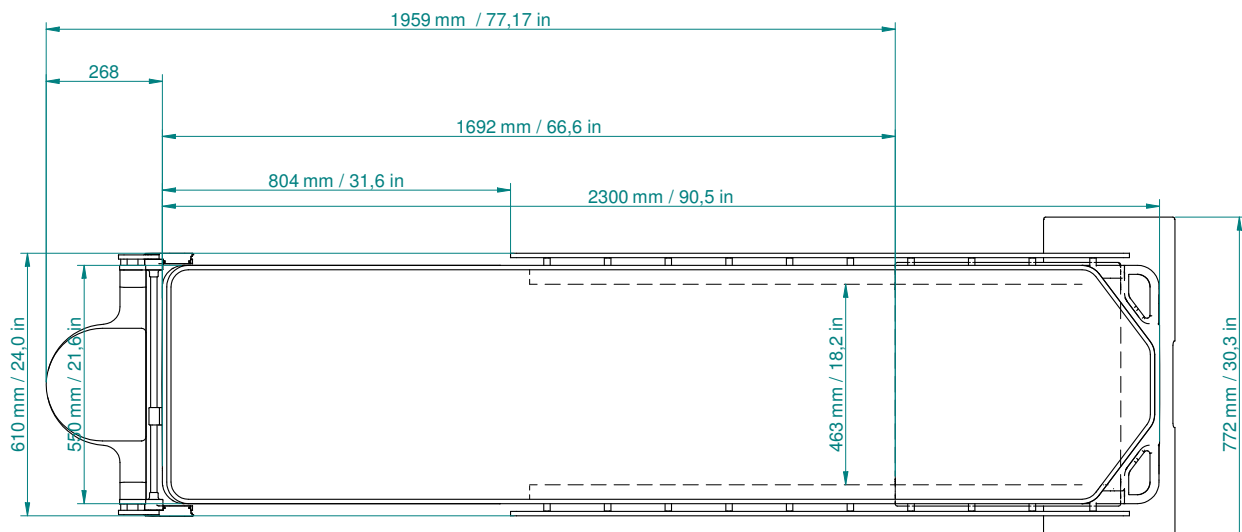

Product specifications of the STILLE imagiQ2

Feature	Description
Patient weight	496 lb / 225 kg and 661 lb / 300 kg (with load upgrade pack)
Table top length , width	L= 90.5" / 230 cm ; W= 21.6" / 55 cm
Translucent length	77.17" / 196 cm with headrest (66.6" / 169.2 cm without headrest)
Translucent table top	Proprietary, patent pending carbon fiber structure. Al 0.4 mm
Float & lock	STILLE True Free Float® and STILLE quick-lock system
Float (longitudinal, lateral, diagonal)	35.4" x 9.9" x 36.8" / 90 cm x 25 cm x 93.4 cm
Operation conditions	100% continuous
AC power: 100-250 Volts	Direct AC power + battery DC power
Length, width, height	L = 92" / 234 cm; W = 30.3" / 77.2 cm; H = 28"- 43" / 71-109 cm
Trendelenburg	±25°
Lateral roll (iso-roll)	±15° (isocentric roll keeps image focus & reduces scattered radiation)
Wheels and steering	4 swiveling back wheels
Side rails	Fixed along 2/3 of the table top
Detachable table top	Yes, future-proof platform for extended use and easy transport.
Transport handle	Yes (Color: dark grey for 225-kg version, silver for 300-kg version)

Customizing the imagiQ2 to fit your needs

Upgrade features	Description
Mobility enhancement pack	When the table is used for frequent intra-clinic transportation, use the mobility enhancement pack for improved mobility. All wheels are designed to swivel and an additional fifth wheel is included to provide optimal mobility. Ref: 535-1711
Load upgrade pack	Maximize the capabilities of the imagiQ2 with increased patient weight capacity up to 661 lb / 300 kg. Ref: 535-1710
Performance pack	For maximum uptime, choose the performance pack. It includes a high-performance long-term battery pack with up to 1 week of durability. Ref: 535-1712